
PROMED

D e n t a l E q u i p m e n t

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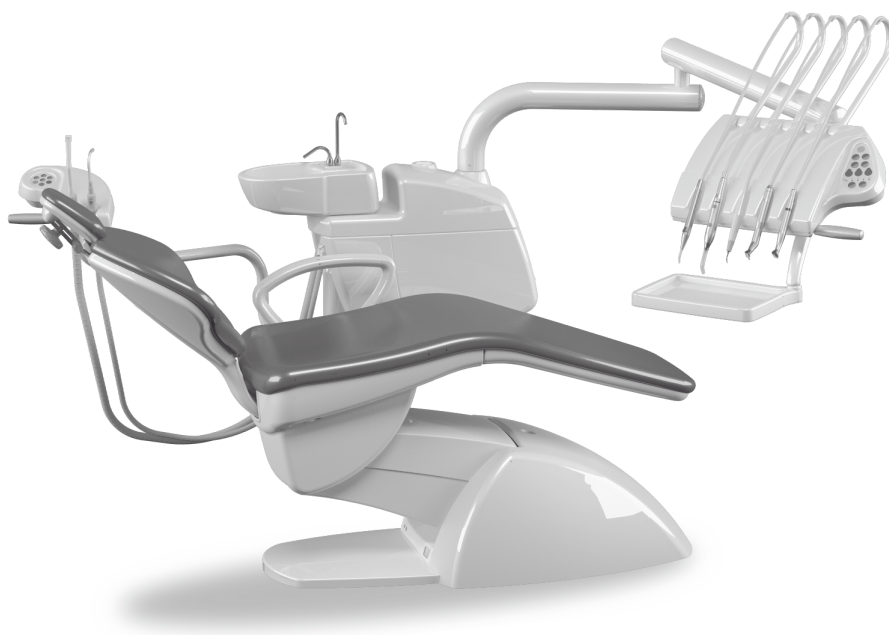Technical data

Trade mark: PROMED
 Model: OASI
 Rated voltage: 230 V
 (~ alternated)
 Rated frequency: 50 Hz
 Rated power: 300 VA
 Product category: Class I Type B
 equipment

Classification according to EN 60 601-1.

5.3 Common equipments (Equipment with cover NOT protected from the water penetration);

5.5 The Equipments (units) are not to be operated in areas with risk of explosion; and closed to any apparatus producing electromagnetic emissions;

5.6 Equipments with continuous operation with intermittent load (=DAB).

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The equipment is conformed to EEC directive 93/42/CEE "Medical Devices"

The manufacturer accepts the responsibility for safety properties of our equipment only if:

- Installation, maintenance, repair and modifications are carried out by a person explicitly authorized by us;
- The general and electric installation of the room is conformed to the related laws;
- The equipment has been correctly set up according to our "operating instructions";
- The equipment has been correctly used according to our "User Manual".

The manufacturer declares that the equipment is a complex dental unit for being used on the man only by a high qualified practical staff in dental practice as job for therapy.

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Attention, see enclosed documentations

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Disposal conforming to the RAEE standard

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Intermittent operation
(30 seconds operation / 300 seconds stop)

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Transportation and storage conditions

When the equipment is packed for shipment and/or storage, it can support the following conditions:

- Environmental temperature from -20°C to +55°C;
- Relative humidity from 5% to 95%; condensation not tolerated.

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- 1 - Auxiliary button
- 2 - Patient chair, up / Patient chair program n°1
- 3 - Patient chair, down / Patient chair program n°2
- 4 - Backrest down / Patient chair program n°3
- 5 - Backrest up / Bowl position, last position
- 6 - Physiological solution to the instruments
- 7 - Predisposition of F.O. / Switch on/off operating lamp
- 8 - Reverse rotation of micromotor / Calling bell
- 9 - Increasing instrument's power / Cup filler
- 10 - Decreasing of instrument's power / Bowl flush

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Assistant's Control panel

- 1 - Patient chair up / Patient chair program n°1
- 2 - Patient chair down / Patient chair program n° 2
- 3 - Backrest down / Patient chair program n° 3
- 4 - Back rest up / Bowl position, last position
- 5 - Switch on/off of operating lamp
- 6 - Bowl flush
- 7 - Cup filler

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Use of the unit

Unit main switch

The unit gets ready for the operating pressing the main switch (position 1). The main switch will be illuminated by a green light.

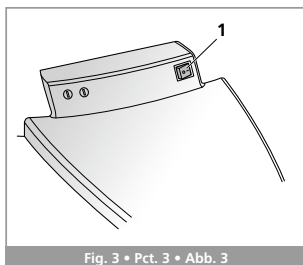
Fig. 3 • Pct. 3 • Abb. 3

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Protection fuse

The dental unit is provided of 2 protection fuses for the power supply. To replace these fuses, you must unscrew the fuse holder (n° 2) and replace the fuse (n° 1). Use only fuses having the same characteristics of the original ones (T 3,15A – 5x20).

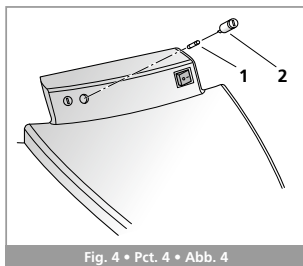
Fig. 4 • Pct. 4 • Abb. 4

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Programming of the patient chair positions

All keys control of the patient chair (table's dentist panel, table's assistant panel and foot control if present) have a double function: if they are pressed (>1 second), they start the driven movements; if they started by an impulse, they started the program.

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The longer pressing key activates the UP movement; with one short pressing key activates program patient chair no. 1.

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The longer pressing key activates the down movement; with one short pressing key activates program patient chair no. 2.

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The longer pressing key activates the inclination (down movement) of the backrest, with one short pressing key activates program patient no. 3.

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The longer pressing key activates the return (up movement) of the backrest. By one short pressing key activates the program that moves to rinsing position. Immediately after reaching this position, automatically the operating lamp is switch off and it's activated the function of the cup filler and then simultaneous the bowl flush function starts. By a second pressing key the patient chair returns to the last work position and then automatically the operating lamp is switched on.

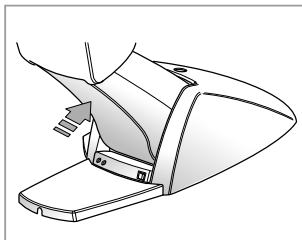
Fig. 5 • Pct. 5 • Abb. 5

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Memorization of programs

For changing one of the other four of the pre-selected functions stored (program no. 1, 2, 3 and LP), bring the patient chair into the desired position manually by pressing the program keys; then press towards the high the lower cover of the arm (Pantograph) of the patient chair and in the same time press the wanted program key and hold depressed until hearing a Beep sound, this is the confirmation of the new storage position instead of the old one.

The new stored program can now be called by pressing the appropriate key.

It is not allowed to move the patient chair when the instruments are taken from holders and are activated by the multifunctional foot control for the keys commands.

DO NOT set up programs when movements are at end stop.

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Cup filling

The cup filling is manual automatic by factory setting; it's activated by pressing key no. 9 of the operator's control panel, only when the instruments are in their holders, or pressing key no. 7 of the assistant's control panel.

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The flushing time can be set adjusting the trimmer on the chair's PCB.

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Bowl Flushing

The flushing time is automatic factory setting; it's activated by pressing key no. 10 of the operator's control panel, only when the instruments are in their holders, or pressing key no. 6 of the assistant's control panel.

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Operating lamp's switch On/Off

The switch on of the operating lamp is activated by pressing key n°7 of the operator's control panel, only when the instruments are in their holders, or pressing key n° 5 of assistant's control panel. A second impulse gives the switch off.

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Key for door opener

Press key no. 8 of the operator's control panel, only when the instruments are in their holders its necessary an external circuit (max 24V-1A) for key for door opener or other services.

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Physiological solution (Optional)

If the bottle for physiological solution is present, you can select by the key board the liquid to the instrument's destination. By the switched on light situated to the upper part of the key board, the distribution of the liquid from the bottle should be made.

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The LED indicator must be switch off during the operating substitution of the bottle.

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Auxiliary push-buttons

The key n° 1 of the dentist's control panel is an auxiliary key, which activation has to be made by an authorized service technician.

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Multifunctional foot control

Pressing the foot lever 1 to the right side, any dynamic selected instrument is started. In case of instruments with a power or adjustable speed (as ultrasonic scaler or micro motor), the lever's opening angle is proportional to the power supplied to the maximum program key.

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Spray management system

Activation of the instrument spray can be configured using two different systems.

STANDARD SYSTEM: Turn lever 1 to turn the instrument on; press lever 2 simultaneously to activate the spray. If lever 2 alone is turned the chip is activated. No command or indication appears on the instruments table.

ALTERNATIVE SYSTEM: The spray must be selected from the instruments table using the "spray tank" button, which has a dual function:

- With the instruments at rest the button controls and, when the LED is on, indicates insertion of sprays from the bottle.
- With the instrument raised the button controls and, when the LED is on, indicates insertion of the spray in the raised instrument.

The spray will be activated when lever 1 on the pedal is turned to start the instrument. If lever 2 is turned the chip is activated. Each instrument stores its last setting. Configuration must be carried out by a qualified technician.

Joystick 3, situated on upper side of the foot control if present, controls the movements of the chair with the same modes as the control panels (see patient's position).

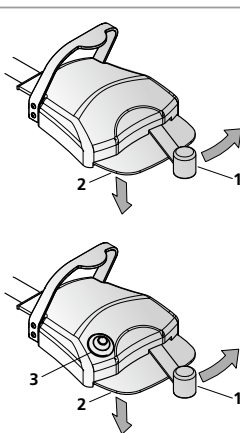
Fig. 6 • Pct. 6 • Abb. 6

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Programmable functions of the instruments

The programming of the function's instruments has been made by the instrument taken from its holder. After returning the same instrument into the holder, all functions are automatically memorized.

When the instrument is taken from its holder, it's possible to set up all desired program functions of the same instrument.

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By the key no. 7 of the dentist's control panel the instrument's F.O. is switched on or switched off; the L.E.D situated at the upper part of the key board shows the predisposition.

The F.O. will be switched on when the foot control will be driven and it will be switched off after 10 seconds from the activation of the same or immediately after returning the instrument into its holder.

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By the key n° 8 of the assistant's control panel, it is switched on or switched off the reverse operation of the micro motor; the switched on LED indicator situated on the upper part of the key board indicates the predisposition of the reverse operation.

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By the key no. 9 and 10 of the dentist's control panel, it's activated the program power or the maximum speed that will be had on the raised instrument, with a foot pedal in activation of the maximum intensity. Such power or speed can be selected on 5 levels: 1% - 25% - 50% - 75% - 100% and will be indicated from the relative LED indicator. Pressing the key + it's possible to increase the speed, while pressing the key - it is possible decrease the speed.

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Instruments

Following instruments are compatible with OASI unit:

- Turbine: all types, if provided with Midwest connection and CE-marked.
- Micromotor: of brand Bien Air models MC3, MC2 and all air micromotors.
- Scaler: of brand Satelec model Suprasson, brand TKD mod. Titanus or brand EMS model Piezo Scaling.
- Curing light: of brand Satelec model Mini LED.

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The pressure regulation of the instruments

Every instrument has an individual regulation of the spray; in the lower part of the dentist's instrument table are placed the regulation taps with the following functions:

- 1 - Screw tap: it regulates the volume of air or the cooling;
- 2 - Tap with Knurl: it regulates the water spray;
- 3 - Screw tap: it regulates the air of chip blower and spray.

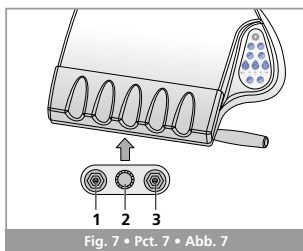
Fig. 7 • Pct. 7 • Abb. 7

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Dentist's syringe and assistant syringe

This instrument does not have the program functions.

- 1 - Control Key for providing water;
- 2 - Control Key for providing air;

Pressing both keys (1 and 2) at the same time it will provide spray. The syringe's cover can be unthread for permitting the sterilization, pressing the realising push button 3.

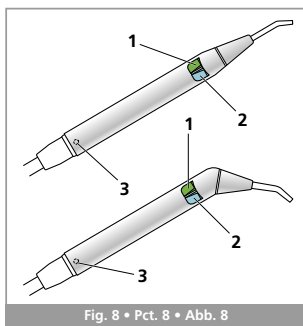
Fig. 8 • Pct. 8 • Abb. 8

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Turbine

When the instrument is taken from its holder, the LED indicates the function's power percentage, automatically it will be appeared as 100% value and it cannot be changed. If the turbine has F.O., it's switched ON (and OFF, resp.) by pressing key board (LED indicator illuminated). The instrument has 3 pressure's regulators situated in the lower part of the back side of the instrument's table.

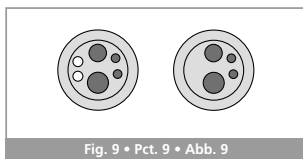
Fig. 9 • Pct. 9 • Abb. 9

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Micromotor

When the instrument is taken from its holder, the same is arranged for the operation. The speed regulation depends on the pressure exercised on the foot control and the key of selection of the spin sense will be active. If the Fibre optic is arranged it will be active the key for switch it on / off. The instrument has the regulator of the spray placed in the lower part of the dentist's instrument table.

The air motors supplied with the unit are equipped by the standard connection type "E", conformed with a directive ISO 3964. It's expected the utilisation of the hand pieces exclusively CE marked.

Electrostatic discharges on the foot control base, can create the stop of the micromotor; in this case it is possible to restore the normal function resetting to zero the foot control (bringing lever 1 of the foot control completely to the left side).

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Ultrasonic Scaler

When the instrument is taken from its holder, the LED indicates the previously memorized power and it can be changed. The instrument has one pressure's regulator for cooling water in the lower part of the doctor's table.

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The oscillation of the ultrasonic scaler to the high frequency can compromise (interfered) the working of the pace-maker.

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Curing Light Mini LED

- 1 - Navigation key for select the working modality
- 2 - Key on/off of working

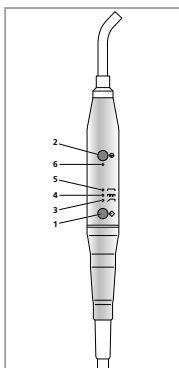
Fig. 10 • Pct. 10 • Abb. 10

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4 - LED indicator "at impulses"

5 - LED indicator "fast"

- 6 - LED indicator working conditions:
Green - normal
Red - Problem of electric alimentation
Red blinking - thermal protection

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Menu “progressive” mode:

select with the navigation key. In “progressive” mode (polymerization similar to menu step with halogen lamp) you will have at disposition:

- Progressive start for 10 seconds ones;
- Full power for 10 seconds one.

Menu “to impulses” mode:

select with the navigation key. This menu corresponds to the full power of the lamp in the mode to impulses, that is an emission of 10 successive flashes with a pause period of 250 Ms in the luminous exposures.

Menu “fast” mode:

after removing the instrument from its holders, the lamp is situated in fast mode. This menu activates the maximum luminous power (according to the terminal part) for 10 seconds.

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For the correct utilisation of the Curing light, it's important to act in the following way:

1. To choose the menu and the lamp is ready for using;
2. Put the terminal part nearest to the surface of the polymerization material to polymerization. It does not need that the terminal part should be in contact with the composite material, as it can be damaged considerably and to decreasing the effect (for the dirt);
3. With one short pressing the key *ON/OFF*, it is switched the polymerization cycle. This action will be confirmed by a “beep” sound.
4. A simple pressure of the key *On/Off* is however possible to interrupt the cycle in any time.
5. If it is not used for 3 consecutive seconds, the curing lamp goes in stand-by and all indicators switched off. It is enough to sufficient press one key for switching it on.

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Optical details

Wavelength: 420-480 nm
Intensity: from 1250 to 2000 mW/cm²
±10% (it depends from the final part)

Tension breaking in the power supply can cause the stopping of the polymerized lamp; in this case it is possible to restore the normal use resetting to zero the foot control (hanging the lever 1 of the foot control to the left completely).

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For the application of others hand pieces, after the first unit's supply, please contact PROMED.

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The headrest can be more or less extracted from the backrest only by pulling or pressing it because it is automatically blocked by a clutch. The double joint for wide articulation can be adjusted by the knob (1): turned in anticlockwise sense, it will be unblocked; in clockwise sense it will be blocked.

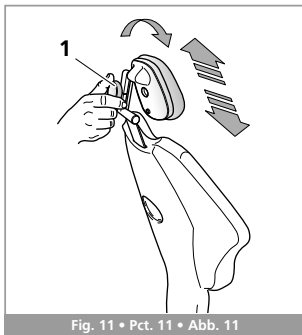
Fig. 11 • Pct. 11 • Abb. 11

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The right armrest can be lowered towards the front side chair to facilitate the patient's access. For doing it turn it forwards and downwards. For restoring it in the normal position make the opposite operation.

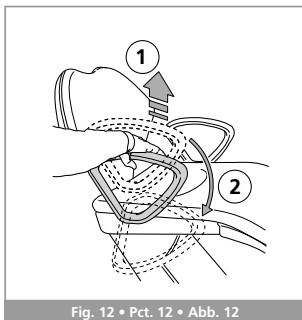
Fig. 12 • Pct. 12 • Abb. 12

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The aspiration should be cleaned daily.

To take it out pull outside the yellow cover and detach the filler making pressure on the snap lock.

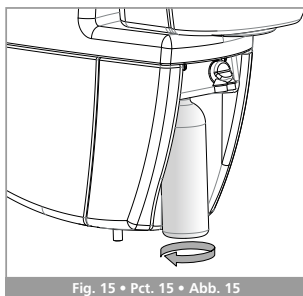
Fig. 15 • Pct. 15 • Abb. 15

Serbatoio alimentazione spray

The alimentation of the spray bottle

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Before taking off the clean bottle, be sure that the LED indicator situated at the upper part of the key board no. 6 of the dentist's instrument table has to be switched off, this guarantees that inside the bottle there is no pressure.

For taking off the bottle, turn it in anticlockwise sense.

To reinstall it repeat the inverse operation.

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For keeping the safety of the device, it is recommended to execute, every one (1) year, the tests of conformity to the electrical safety.

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Cleaning and disinfection of the unit

Surface cleaning

It is advised to clean up the external surfaces of the unit with hygienic liquid for a plastic material.

The instruments tubes must be washed with water and neutral soap, then well dry.

The upholstery of the patient chair and the stools have to be cleaned up with water and neutral soap.

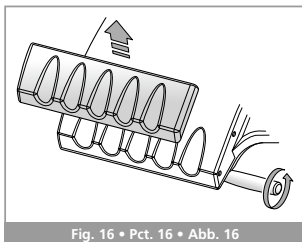
Fig. 16 • Pct. 16 • Abb. 16

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Handles and protection in silicone

The handles of the dentist's and assistant instrument table are taken off turning in clockwise sense and are cleaned by a disinfectant liquid.

The protection in silicone of the instruments table can be taken off for a better cleaning and disinfection using a disinfectant liquid.

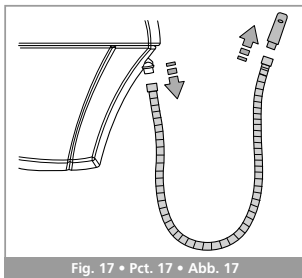
Fig. 17 • Pct. 17 • Abb. 17

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Cannulas and aspiration tubes

The connectors and the aspiration tubes can be taken off for cleaning up and disinfecting, from the relative connectors.

The cleaning of the cannulas must be done by disinfectant liquid.

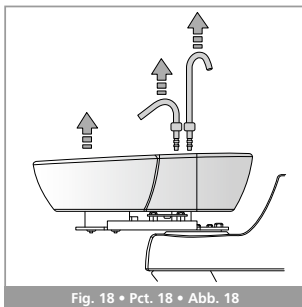
Fig. 18 • Pct. 18 • Abb. 18

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Cuspidor bowl

The cuspidor and the nozzles of the cup filler and washing can be removed for cleaning and disinfection pulling them. The cleaning of the small items must be done by a disinfectant liquid.

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Product Identification

For any communication with the manufacturer or after sales centres, it must be always mentioned the equipment's model (mod.) and the serial number printed on the matriculation label situated on the arm of the installation.

Technical notes

The wiring diagrams of the equipment are available in the qualified technical central services. The technical assistance on the equipment must be done exclusively by a high qualified technician.

Disposing of the device

Definition of RAEE

The symbol "REEE" means Refusal of Electrical and Electronically Equipment; it has to be considered as RAEE all electrical and electronically equipments for domestic or professional use that at the end of their life must be "scraped".

Disposing modality of the equipment

At the end of the usefull life, expected by the manufacturer in 10 years from the installation date, the equipment must be made unusable following this instructions:

- Untie the electrical, hydraulic and pneumatic connections of the equipment;
- Cut all electrical connections which go from the base of the equipment to different peripherals;
- Cut the cables of the secondary exit of the transformer as they cannot be restored.

When the equipment is useless, it has to be packed by recyclable material and as to be decided from the customer or the retailer, it has to be sent to private or public rubbish dumps equipped for the disposing of the REEE. The customer or the retailer must take trace of the made "Disposing" of the equipment through the suitable structures, keeping the transport documents with precise motive's indication.

Attention!! The symbol placed on the equipment, does not permit the disposing of the same with different ways from what we explained above.

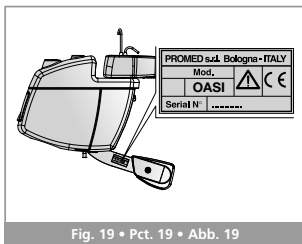
Fig. 19 • Pct. 19 • Abb. 19

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Enclosure I: Data according to EN 60601-1-2 (2003)

OASI Dental Unit has been manufactured for use in an environment like the one cited below. The customer or user of OASI should ensure that it is used in the correct environment.

Transmission Tests	Conformance	Electromagnetic environment
RF transmission	Group 1	OASI use RF energy for its internal operation. Its transmission is very low and it is highly probable that it doesn't disturb neighbouring electronic devices.
RF transmission	Class b	OASI is for use in all premises including residential ones and all premises that are, directly connected to a public low voltage power supply covering residential buildings.
Harmonic transmission	Class A	
Transmission of voltage Fluctuations or Flicker	Conforms	

Immunity test	EC 60601 test level	Conformance level	Electromagnetic environment
Electrostatic discharge	± 6 kV contact discharge ± 8 kV atmospheric discharge	± 6 kV contact discharge ± 8 kV atmospheric discharge	Floors made of wood or concrete or have ceramic tiles. If they are covered with synthetic material, the relative humidity must be at least 30%
Fast transient electrical disturbances / lines	± 2 kV for power lines ± 1 kV for input and output	± 2 kV for power lines ± 1 kV for input and output	The quality of the supply should correspond to that of a typical business or hospital environment
Surges	± 2 kV for power lines ± 1 kV for input and output	± 2 kV for power lines ± 1 kV for input and output	The quality of the supply should correspond to that of a typical business or residential environment
Voltage interruptions, short-term interruptions and fluctuations of the supply voltage	<5% Ut (95% interruption) for 0,5 period 40% Ut (60% interruption) for 5 periods 70% Ut (30% interruption) for 25 periods <5% Ut (95% interruption) for 5 seconds	<5% Ut (95% interruption) for 0,5 period 40% Ut (60% interruption) for 5 periods 70% Ut (30% interruption) for 25 periods <5% Ut (95% interruption) for 5 seconds	The quality of the supply should correspond to that of a typical business or residential environment. If the user of OASI asks continued operation even if the power supply is interrupted, it is recommended to supply OASI Dental Unit from an uninterrupted power supply.
Magnetic field with a supply (50Hz)	3 A/m	3 A/m	Magnetic fields at the mains frequency should correspond to typical values in a business and residential environment.

Operating instruction

Immunity test	EC 60601 test level	Conformance level	Electromagnetic environment
Conducted HF disturbances	3 Veff	3 Veff	$d = 1,17\sqrt{P}$ with P as the maximum rated power of the transmitter in Watts (W) according to the transmitter manufacturer and d as the recommended safe distance in meters
Radiated HF disturbances	3 V/m	3 V/m	$d = 1,17\sqrt{P}$ for 80 Mhz to 800 Mhz $d = 2,33\sqrt{P}$ for 800 Mhz to 2,5 Ghz with P as the maximum rated power of the transmitter in Watts (W) according to the transmitter manufacturer and d as the recommended safe distance in meters

OASI is intended for use in an electromagnetic environment in which HF disturbances are controlled. The user of OASI can help prevent electromagnetic disturbances by maintaining the minimum distance between portable and mobile HF telecommunications devices (transmitters) and the OASI depending on the output of the communication device as indicated.

Rated power of the transmitter in W	Safe distance depending on the transmission frequency:			
	150 KHz to 80 Mhz	80 Mhz to 800 Mhz	800 Mhz to 1 Ghz	1 Ghz to 2,5 Ghz
0,01	0,12	0,04	0,07	0,23
0,1	0,37	0,11	0,22	0,74
1	1,17	0,35	0,70	2,33
10	3,69	1,11	2,21	7,38
100	11,67	3,5	7,00	23,33

Maximum length of power cords:

- Chair supply: 1.2 mt AC
- Foot control: 1.5 mt

WARNING: The use of accessories, transducers and cables other than those specified, except those sold by the manufacturer as spare parts, may cause increased emissions or decreased immunity of the device.

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